



Safety Instructions Please read this manu it properly. **▲** Warning Attention: • It is forbidden to use this product in the head, face, elbow, knee and other joints, spine and other bone parts with less muscle coverage. cervical vertebrae of the neck, so as to avoid damaging the carotic sinus, resulting in fainting and cardiac arrest.

- It is forbidden to use this product in the armpit, the medial side o upper arm, the medial side of clavicle, the medial side of elbow joint, the popliteal fossa at the back of knee, abdomen, kidney area, genitalia
- Patients with kidney stones, ureteral stones, bladder stones, urinar tract lithiasis, urinary tract infection and other urinary diseases are prohibited to use this product below the abdomen. • It is forbidden to use this product in the body area with bruise, pa surface wound, local fever, swelling discomfort, acute inflammatic muscle spasm, muscle laceration, acute sprain, contusion and strain • In case of fever, nausea, vomiting, muscle pain and weakness after strenuous exercise, it is forbidden to use this product and you shoul

and other vulnerable parts with dense nerves and blood vessels.

seek medical advice immediately. • Patients with the following diseases are prohibited to use this produce cardiovascular and cerebrovascular diseases, blood diseas coagulation dysfunction, thrombocytopenic purpura, neck plaqu varicose veins, venous thrombosis, cancer, tumor, mental diseases. • Patients with fracture or internal nail are prohibited to use this produ • It is forbidden to use this product in precise operation or driving. • Please do not use this product continuously for a long time in order prevent muscle damage caused by excessive massage.

Safety Instructions

and fatigue, relieve muscle pain and discomfort, and improve you professional medical equipment. It has no effect on weight loss, ca not cure any disease, and can not solve medical problems. There or improper operations.

Dos and Don'ts • The product has been charged in advance with about 60% to 80% o

the power before leaving the factory, and can be used without beir charged after opening. · When the machine is running, do not put fingers, hair or other par of the body close to the interface between the massage head and th machine so as to avoid injury. of massage head and the charging interface, otherwise it may cau scald, electric shock or fire due to short circuit. Please do not soak the product in water or other liquids. Please do not use this product in humid environment.

Please use and charge it at 0 ~ 40 °C.

affect battery performance. · Please do not repair, disassemble or modify the product by yourse Please do not use non-original charger for charging, and only use the very original charger provided by our company. Please do not use this product when it is being charged. After finishing charging, please take away the charger from th • If you have not used this product for a long time, please charge it in advance before you reuse it someday.

Keep away from water, fire and corrosive environment, so as no

Safety Instructions

Insuitable people spital before use, and only use with permission.

People over 60 years old. People with limb loss or/and high paraplegia. People with pacemaker and stent.

eople who took some surgery within 12 months. eople with orthopedic diseases such as abnormal spine, joints an bones, osteoporosis, cervical spondylosis, lumbar disc herniatio disease or/and age. neuropathy or retinal injury. Patients with kidney stones, ureteral stones, bladder stones, urinary

n case of abnormality and failure:

immediately and contact the after-sales service center.

Product introduction

nder age in oder to avoid accidents. If any of the following conditio e met, please consult the surgeon or attending doctor of a 3.

tract lithiasis, urinary tract infection and other urinary diseases.


Avoid beating, shaking or hitting the product. Strong impact may caus damage to parts of the machine and lead to abnormal work. In case

nperature of the machine to rise and damage the machine. Please do not use it continuously for a long time. In case of fever, smoking, excessive buzzer of the machine, please p using it immediately and contact the after-sales service center high temperature, deformation, expansion and battery leakage



tight muscles.



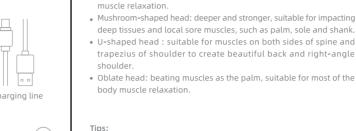

To relieve the back and shoulder muscl

by sedentary and bad posture.

Silicone massage heads

Daily massage to relax the body.

Product Function Introduction

To relieve the delayed muscle soreness after exercise and relax t

· Round Bobo head: Soft and elastic, suitable for most of the boo

are schematic diagrams for reference only. Due to the updating an development of the product, there may be some differences between the actual product and the schematic diagram. Please refer to the actual

Product Function Introduction

uscle injury caused by prolonged massage, the gear indicator light nt remains in the state. Click the power key to enter the first gear hin 5 seconds to time it again. If there is no any operation, the achine will shut down automatically.

cked-Rotor protection: sts for more than 3 seconds, it will enter standby mode. Click power y to reactivate in 10 seconds. If there is no any operation, it will shu

Maintenance



Please install and remove the massage head in the state of power o and ensure that the head is installed firmly before starting the machin When using it, the massage head cannot be adjusted or loaded so a

• If you have not used this product for a long time, please charge it

 If high temperature, deformation, expansion and battery leaka occur in the process of using the machine, please stop using immediately and contact the after-sales service center.

advance before reusing it.

Installation process

The Fascia Massage Gun has its own massage head (flat head), which can be used independently or combined according to the demand. The installation effect is as follows.



it off and charge it in time.

Usage

1. Indicator light

Gear indicator (gear shift) Gear level Gear indicator (breathing light) First gear Blue light breathing 1500 Arousing muscle

Second gear Green light breathing 2400 Massage to relax

Third gear Red light breathing 3200 Deep massage

Tester: 08 Test result: Qualified

Shenzhen YUNMAI Technology Co., Ltd.

Standby gear (shut down after 2s) 0 Power indicator (power indicator in use) • The green light is on in a stable state, which means full power. The green light flashes slowly, which means medium power. • The green light flashes guickly, which means low power. Please

Fourth gear Yellow light breathing 650

ong press the power key to turn it on or off. After power on, single click the power key to change gears. It cou ycle in the order of one, two, three and four gears in turn. is recommended to run it for 5 seconds after starting the machi ore using it. It is forbidden to use it when it is being charged!

lassage Gun . Insert the charging line into the charging port and connect the power supply to start charging. Power indicator (charging indicator) • It is being charged if the red light is on. • The charging is completed with full power if the green light is or ▲ Attention


Do not use it when it is being charged!

to your needs. And it is recommended to try it from the low gear level Gently press the massage head on the muscle part to be relaxed, and area, apply appropriate pressure at the tenderness point, and stay for to use excessive force when applying pressure so as to avoid damaging in the same part for more than 60s. he massage time appropriately. and the continuous use time should not exceed 30 minutes. • Please operate the product in strict accordance with the massage position and operation methods written in this manual.

Massage Gun , then you can change to the appropriate gear according

Usage

4. How to use

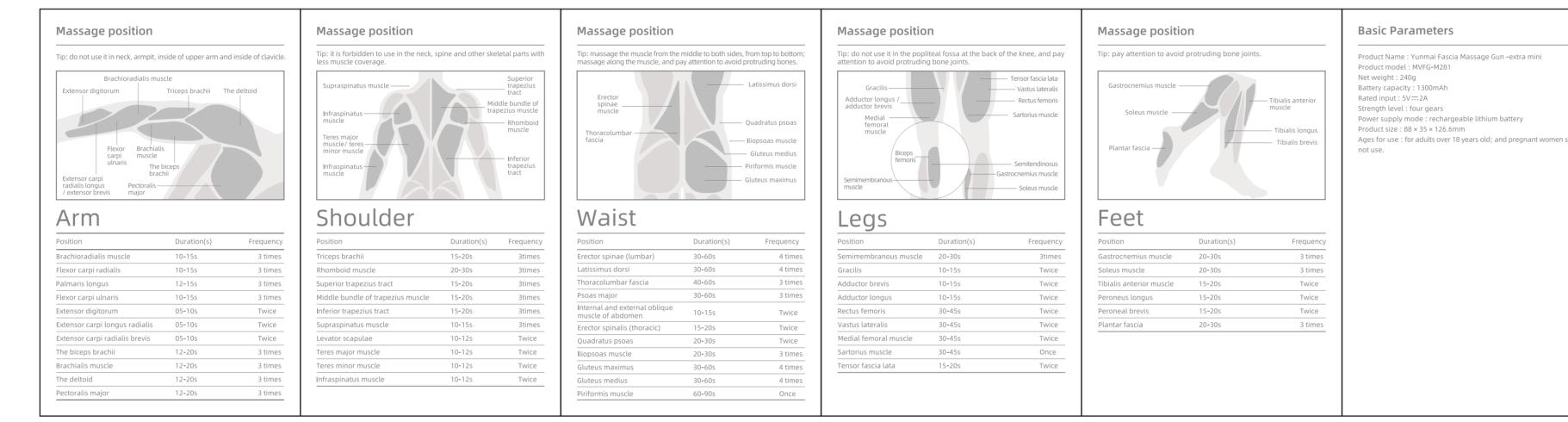



1. Cleaning • Please wipe the machine body and massage heads with clean and slightly wet soft cloth and wet towel regularly; then dry them with dry cloth and paper towel, and finally put them in the recycling bag or Certificate • When cleaning, make sure that the power switch is turned off and th Thank you for purchasing our product! charging plug is removed.

• Please do not use corrosive solution to wipe the machine bod Product name: Yunmai Fascia Massage Gun extra min discoloration. Product model: MVFG-M281 • Please store the product in a cool and dry place. Please do not store the machine body in high temperature environments Through strict examination and quality inspection of away from humid, water, fire and corrosive environment, so as r the above product by our professional quality inspectors to affect the battery performance. we hereby confirm that the product passed the · When it needs to be stored or carried out, the Fascia Massage G inspection standards of Shenzhen Yunmai Technology can be stored in the attached portable storage bag.

• Please do not touch, squeeze, drop, throw or/and scratch the product etc. During transportation, the product should be fixed steadily a put up gently. We should prevent the product from extrusion, viole

备注: 首页和尾页背面不印刷。







